



A GUIDE TO THE PROPERTIES AND USES OF

Detergents

IN BIOLOGY AND BIOCHEMISTRY

*The
Detergent
Standard
for over
50 Years*



DETERGENTS

**A guide to the properties
and uses of detergents
in biological systems**

By

Srirama M. Bhairi, Ph.D.

A Word To Our Customers

CALBIOCHEM® has been the leading supplier of a variety of quality detergents to a large number of researchers all over the world for almost 50 years. During this period, we have received a number of inquiries on the use of detergents, definitions, relevance of critical micelle concentration (CMC), cloud point, hydrophilic number, and how to select the most appropriate detergent. As a service to the research community, CALBIOCHEM® is providing this guide to the use of detergents in biological systems. The background information and the selected bibliography provided here will hopefully serve the needs of the first time users of detergents as well as those of experienced investigators.

We have also included a chapter on a unique series of compounds known as Non-Detergent Sulfobetaines (NDSBs). As evident from the name, these compounds are not detergents and they do not form micelles. Structurally, NDSBs have hydrophilic groups similar to those found in zwitterionic detergents; however, they possess a much shorter hydrophobic chain. They have been reported to improve the yield of membrane proteins when used in conjunction with the traditional detergents and prevent aggregation during renaturation of chemically or thermally denatured proteins.

The discussion provided in this booklet is by no means complete. However, we hope it will help in the understanding of general principles involved in the use of detergents.

Margaret Dentlinger
Director of Marketing

Table of Contents

| | |
|--|-----------|
| Hydrophobic Interactions | 4 |
| What are Detergents? | 6 |
| Biological Membranes | 7 |
| How do Detergents Solubilize Membrane Proteins? | 8 |
| Classification of Detergents | 10 |
| General Detergent Properties | 12 |
| Removal of Unbound Detergents | 15 |
| Guidelines For Choosing a Detergent | 17 |
| Appendix | 19 |
| Non-Detergent Sulfobetaines | 19 |
| CALBIOSORB™ Adsorbent | 20 |
| Structure and Classification of Detergents | 24 |
| Selected Bibliography | 26 |
| CALBIOCHEM® Detergents | 35 |

Hydrophobic Interactions

Water forms a highly ordered network of intermolecular hydrogen bonds (Figure 1). It is the strength of all the hydrogen bonds combined that imparts the liquid properties to water. Polar or hydrophilic substances dissolve in water because they form hydrogen bonds and electrostatic interactions with water molecules. Non-polar or hydrophobic substances, on the other hand, are unable to form such interactions, and consequently, are immiscible with water. Addition of nonpolar substances to water disrupts intermolecular hydrogen bonding of water molecules and creates a cavity which is devoid of the water molecules. At the surface of the cavity, water molecules rearrange in an orderly manner (Figure 2). This results in a thermodynamically unfavorable decrease in entropy. To compensate for the loss of entropy, water molecules force the hydrophobic molecules to cluster and thus occupy the minimum space. This phenomenon is known as the **hydrophobic effect** and the “forces” between hydrophobic regions are called **hydrophobic interactions**.

Hydrophobic interactions play a major role in defining the native tertiary structure of proteins. Proteins consist of polar and non-polar amino acids. In water-soluble proteins, hydrophobic domains rich in non-polar amino acids are folded in together and thus are shielded from the aqueous environment. In membrane proteins, some hydrophobic regions that otherwise would be exposed to the aqueous environment are surrounded by lipids.

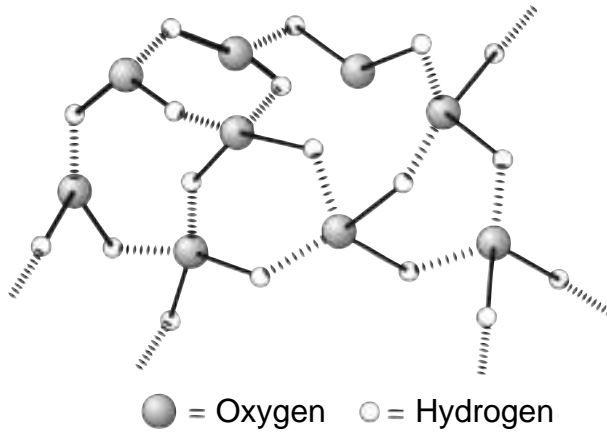


Figure 1: Inter-molecular hydrogen bonding in water.

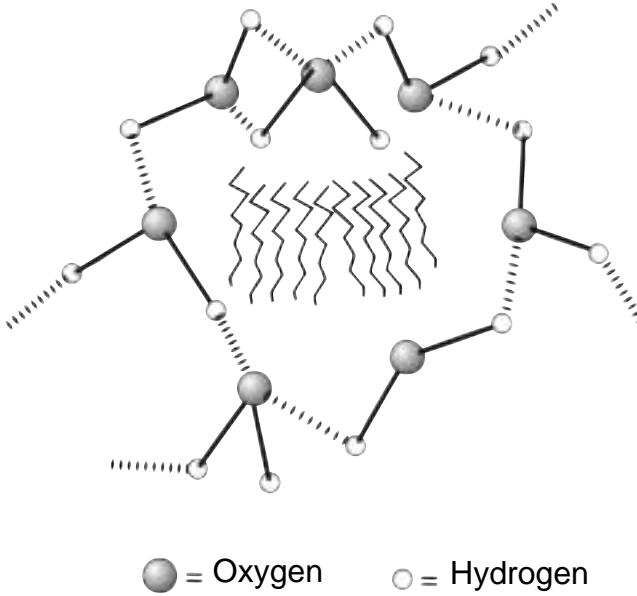


Figure 2: Clustering of hydrocarbon molecules in water.

What are Detergents?

Detergents are amphipathic molecules that contain both polar and hydrophobic groups. These molecules contain a polar group (head) at the end of a long hydrophobic carbon chain (tail). In contrast to purely polar or non-polar molecules, amphipathic molecules exhibit unique properties in water. Their polar group forms hydrogen bonds with water molecules, while the hydrocarbon chains aggregate due to hydrophobic interactions. These properties allow detergents to be soluble in water. In aqueous solutions, they form organized spherical structures called **micelles** (Figure 3), each of which contain several detergent molecules. Because of their amphipathic nature, detergents are able to solubilize hydrophobic compounds in water. Incidentally, one of the methods used to determine the CMC (see below) relies on the ability of detergents to solubilize a hydrophobic dye. Detergents are also known as surfactants because they decrease the surface tension of water.

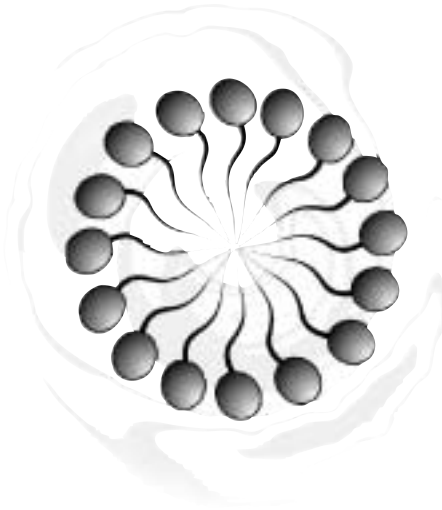


Figure 3: A detergent-micelle in water.

Biological Membranes

Biological membranes are composed of phospholipids and proteins where phospholipids can be viewed as biological detergents. The majority of the lipids that make up the membrane contain two hydrophobic groups connected to a polar head. This molecular architecture allows lipids to form structures called **lipid bilayers**, in which the hydrophobic chains face each other while the polar head groups are outside facing the aqueous milieu (Figure 4). Proteins and lipids, like cholesterol, are embedded in this bilayer. This bilayer model for membranes was first proposed by Singer and Nicolson in 1972 and is known as the **fluid mosaic model** (Figure 5). The embedded proteins are held in the membrane by hydrophobic interactions between the hydrocarbon chains of the lipids and the hydrophobic domains of the proteins. These membrane proteins, known as **integral membrane proteins**, are insoluble in water but are soluble in detergent solutions.

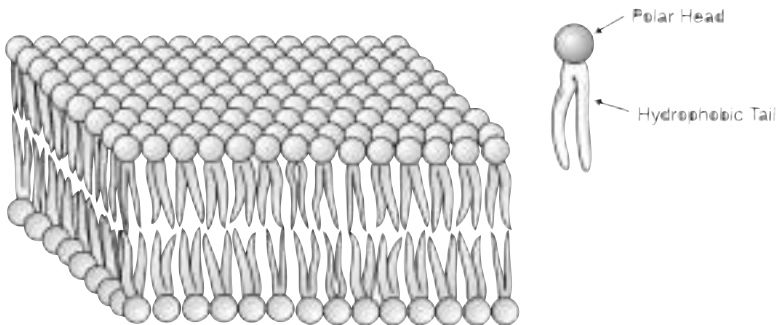


Figure 4: A phospholipid bilayer.

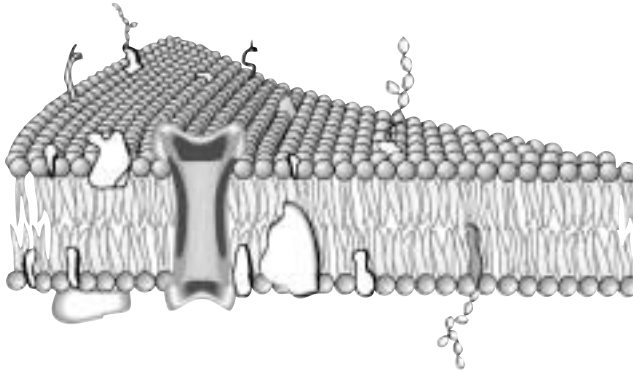


Figure. 5: Fluid-mosaic model of a biological membrane.

How Do Detergents Solubilize Membrane Proteins?

Detergents solubilize membrane proteins by mimicking the lipid-bilayer environment. Micelles formed by detergents are analogous to the bilayers of the biological membranes. Proteins incorporate into these micelles via hydrophobic interactions. Hydrophobic regions of membrane proteins, normally embedded in the membrane lipid bilayer, are now surrounded by a layer of detergent molecules and the hydrophilic portions are exposed to the aqueous medium. This keeps the membrane proteins in solution. Complete removal of detergent could result in aggregation due to the clustering of hydrophobic regions and, hence, may cause precipitation of membrane proteins.

Although, phospholipids can be used as detergents in simulating the bilayer environment, they form large structures, called vesicles, which are not easily amenable for isolation and characterization of membrane proteins. Lyso-phospholipids form micelles that are similar in size to those formed by many detergents. However, they are too expensive to be of general use in everyday protein biochemistry. Hence, the use of synthetic detergents is highly preferred for the isolation of membrane proteins.

Dissolution of membranes by detergents can be divided into different stages (Figure 6). At low concentrations, detergents bind to the membrane by partitioning into the lipid bilayer. At higher concentrations, when the bilayers are saturated with detergents, the membranes disintegrate to form mixed micelles with the detergent molecules. In the detergent-protein mixed micelles, hydrophobic regions of the membrane proteins are surrounded by the hydrophobic chains of micelles. In the final

stages, solubilization of the membranes leads to the formation of mixed micelles consisting of lipids and detergents and detergent micelles containing proteins (usually one protein molecule per micelle). For example, solubilization of a membrane containing rhodopsin by digitonin leads to complexes containing one rhodopsin molecule per micelle consisting of 180 digitonin molecules. Other combinations of micelles containing lipids and detergents and lipid-protein-detergent molecules are possible at intermediate concentrations of detergent. Micelles containing protein-detergent molecules can be separated from other micelles based on their charge, size, or density.

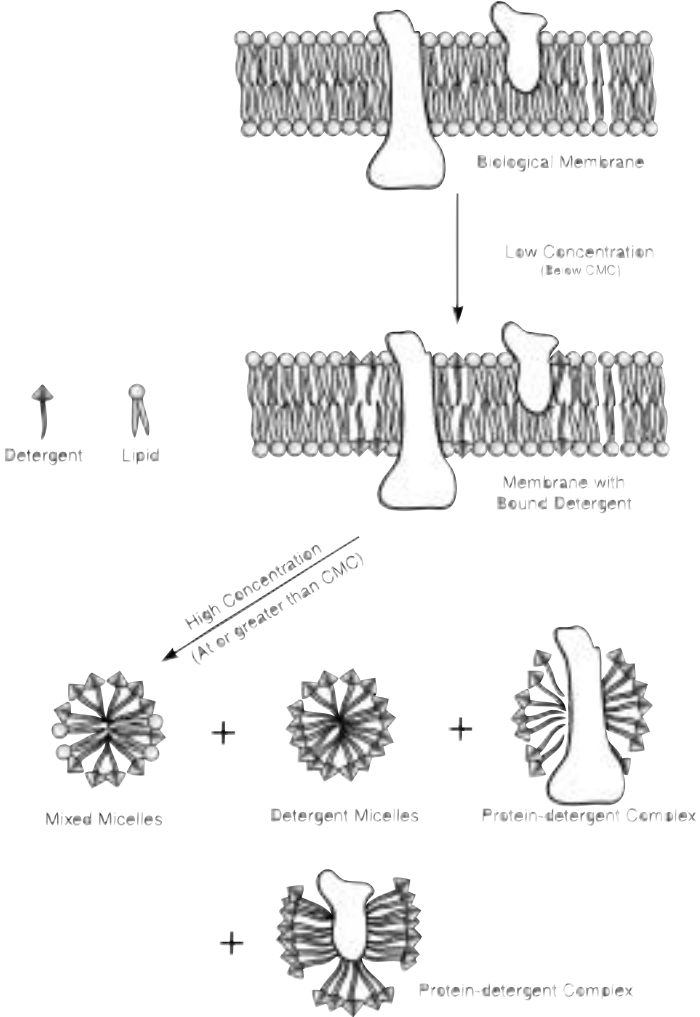


Figure 6: Stages in the dissolution of a biological membrane with detergents.

Classification of Detergents

A large number of detergents with various combinations of hydrophobic and hydrophilic groups are now commercially available. Based on the nature of the hydrophilic head group, they can be broadly classified as **ionic**, **non-ionic**, and **zwitterionic** detergents.

Ionic Detergents

Ionic detergents contain a head group with a net charge. They can be either negatively (anionic) or positively charged (cationic). For example, sodium dodecyl sulfate (SDS), which contains the negatively charged sulfate group, is an anionic detergent while cetyl trimethyl-ammonium bromide (CTAB), which carries the positively charged trimethylammonium group, is a cationic detergent. Furthermore, the ionic detergents either contain a hydrocarbon (alkyl) straight chain as in SDS and CTAB, or a more complicated rigid steroidal structure as in sodium deoxycholate (see bile acid salts). There is a repulsion between the similarly charged polar groups of detergent molecules in a micelle. Therefore, the size of the micelle is determined by the combined effect of hydrophobic attraction of the side chains and the repulsive forces of the ionic groups. Consequently, neutralizing the charge on the head group with increasing concentrations of a counter ion leads to a larger micellar size. Micellar size also increases with the increase in alkyl chain length.

Bile Acid Salts

Bile acid salts are anionic detergents containing a rigid steroidal hydrophobic group (e.g., sodium salts of cholic acid and deoxycholic acid). In addition to the anionic carboxyl group at the end of the short alkyl chain they also carry hydroxyl groups on the steroid structure. Thus, there is no well-defined polar head group. Instead, the bean shaped molecule has a polar and an apolar face.

Bile acid salts form small aggregates. They can be conjugated to taurine or glycine at the end of the carboxyl group. Unlike spherical micelles formed by alkyl ionic detergents, the micelles formed by bile acid salts are kidney shaped due to their rigid structure. As for ionic detergents, their micellar size is influenced by the concentration of the counter ion. Due to the low pK_a (5 - 6) of the unconjugated bile salt, and low solubility of bile acids, their utility is limited to the alkaline pH range. On the other hand, the pK_a of conjugated bile acid salts is much lower, hence, they can be used over a broad pH range. Dihydroxy bile acid salts and deoxycholate are more effective than trihydroxy bile acid

salts in membrane solubilization and in dissociation of protein-protein interactions. Trihydroxy bile acid salts are milder and are better suited for removal by dialysis.

CALBIOCHEM® offers over 20 anionic detergents. Among them you will find BATC (Cat. No. 196950) and TOPPS (*tert*-octyl phenyl propanesulfonic acid) (Cat. No. 615000). BATC (4'-amino-7-benzamido-taurocholic acid), a synthetic derivative of taurocholic acid, is especially suited for solubilization of glycosyl-phosphatidyl inositol (GPI) anchored membrane proteins. TOPPS has the same aromatic hydrophobic tail as TRITON® X-100; however, instead of non-ionic polyoxyethylene groups, it contains an ionic sulfonate group as the hydrophilic head group. This detergent has been used in the renaturation and refolding of chemically and thermally denatured carbonic anhydrase B.

Non-ionic Detergents

Non-ionic detergents contain uncharged, hydrophilic head groups that consist of either polyoxyethylene moieties as in BRIJ® and TRITON® or glycosidic groups as in octyl glucoside and dodecyl maltoside. In general, non-ionic detergents are better suited for breaking lipid-lipid and lipid-protein interactions than protein-protein interactions. Hence, they are considered non-denaturant and are widely used in the isolation of membrane proteins in their biologically active form. Unlike ionic detergents, salts have minimal effect on the micellar size of the non-ionic detergents.

Detergents with polyoxyethylene head groups, may contain alkylpolyethylene ethers with the general formula $C_nH_{2n+1}(OCH_2CH_2)_xOH$, or a phenyl ring between the alkyl chain and the ether group. TRITON® X-100 and NP-40 belong to the latter class (see Table I). Polyoxyethylene chains form random coils and are consequently removed farther from the hydrophobic core of the micelles. Detergents with shorter polyoxyethylene chains form aggregates and viscous solutions in water at room temperature, whereas those with longer chains do not aggregate. It should be noted that detergents containing aromatic rings absorb in the ultraviolet region. They may interfere with spectrophotometric monitoring of proteins at 280 nm. CALBIOCHEM® offers hydrogenated versions of these detergents in which the aromatic rings are reduced and these detergents exhibit relatively low absorption at 280 nm.

Alkyl glycosides have become more popular as nonionic detergents in the isolation of membrane proteins for several reasons. First, they are homogeneous with respect to their composition and structure. Second, several variations of alkyl glycosides containing different combinations of the hydrocarbon chain (cyclic or straight chain) and the polar sugar

group can be easily synthesized in pure forms. Third, subtle differences in the physicochemical properties of alkyl glycosides bearing various alkyl chains, attached to either to a glucose, maltose, or a sucrose head group, can be exploited for selective solubilization of membrane proteins.

Zwitterionic Detergents

Zwitterionic detergents are unique in that they offer the combined properties of ionic and non-ionic detergents. Like non-ionic detergents the zwittergents, including CHAPS and the ZWITTERGENT® 3-X-series, do not possess a net charge, they lack conductivity and electrophoretic mobility, and do not bind to ion-exchange resins. However, like ionic detergents, they are efficient at breaking protein-protein interactions. Zwittergents such as CHAPS are less denaturing than the ZWITTERGENT® 3-X series, possibly owing to their rigid steroid ring structure.

General Properties of Detergents

Critical Micelle Concentration (CMC)

The CMC can be defined as the lowest concentration above which monomers cluster to form micelles. Alternatively, it is the maximum attainable chemical potential (concentration) of the monomer. In reality, micellization occurs over a narrow concentration range rather than at a particular concentration. The CMC decreases with the length of the alkyl chain and increases with the introduction of double bonds and other branched points such as would occur in bile acid salts. Additives, such as urea, that break up water structure also increase the CMC. In ionic detergents, the CMC is reduced by increasing the concentration of counter ions. From a practical point of view, a high CMC is desirable when dialysis is used for the removal of the detergent.

Three of the most popular methods used to determine CMC are surface tension, light scattering, and dye solubilization. Surface tension decreases with the detergent concentration and reaches a minimum around the CMC value. Light scattering as well as the solubility of a hydrophobic dye increase with detergent concentration. The point of inflection on a graph obtained by plotting any of the three parameters vs the detergent concentration corresponds to the CMC of the detergent (Figure 7).

Given the CMC, the concentration of the detergent and the aggregation number (see below), it is possible to calculate the concentration of

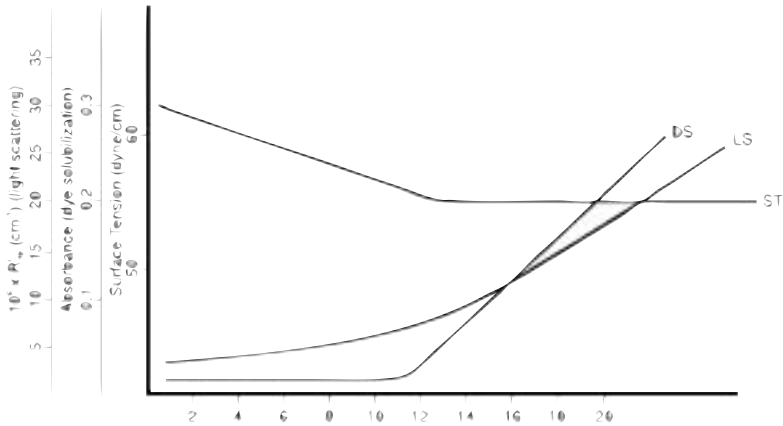


Figure 7: Representative results for determining the CMC of a surfactant by various methods.

Note: DS = dye solubilization; LS = light scattering; ST= surface tension.

micelles in moles per liter using the following formula:

$$[\text{micelles}] = (C_s - \text{CMC}) \div N$$

where C_s is the bulk molar concentration of detergent and N is the mean aggregation number. For example, a solution containing 35 mM of CHAPSO (M.W. = 630.9) in PBS buffer will have $[(35 - 8) \div 11]$ or 2.45 mM of micelles.

Kraft Point

At very low temperatures, detergents remain mainly in an insoluble crystalline state and are in equilibrium with small amounts of dissolved monomer. As the temperature increases, more and more of the monomeric detergent goes into solution until the concentration of the detergent reaches the CMC. At this point it exists predominantly in the micellar form. The temperature at which the monomer reaches the CMC concentration is called *critical micellar temperature* (CMT). The temperature at which all the three phases - crystalline, micellar and monomeric - exist in equilibrium is called the *Kraft Point* (Figure 8). At this temperature the detergent solution turns clear and the concentration of the detergent reaches its CMC value. For most detergents, the Kraft point is identical to the CMT.

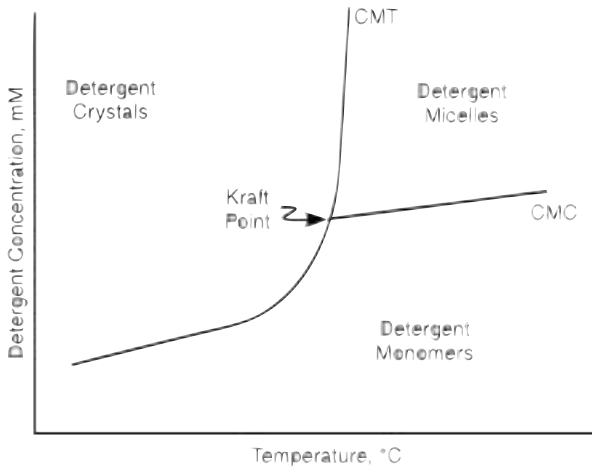


Figure 8: Temperature-composition phase diagram for detergent solutions.

Cloud Point

At a particular temperature above the CMT, non-ionic detergents become cloudy and undergo phase separation to yield a detergent-rich layer and an aqueous layer. This temperature is called the cloud point. Phase separation presumably occurs due to a decrease in hydration of the head group. For example, the *Cloud Point* of TRITON® X-100 is 64°C whereas that for TRITON® X-114 is around 22°C. Hence, TRITON® X-114 solutions are maintained cold. This property can be used to a particular advantage. Membranes can be at first solubilized at 0°C and the solution can be warmed to about 30°C to effect the phase separation. This allows partition of integral membrane proteins into the detergent rich phase which can be later separated by centrifugation.

Aggregation Number

This is the number of monomeric detergent molecules contained in a single micelle. It can be obtained by dividing the micellar molecular weight by the monomeric molecular weight. The molecular weight of micelles can be obtained from various techniques including gel filtration, light scattering, sedimentation equilibrium, and small-angle X-ray scattering. The micelles formed by bile acid salts tend to have low aggregation numbers while those formed by TRITON® have high aggregation numbers. Like micellar size, the aggregation number is also influenced by the ionic strength.

$$\frac{\text{micellar molecular weight}}{\text{monomeric molecular weight}} = \text{aggregation number}$$

Hydrophile-Lipophile Balance (HLB)

This is a measure of hydrophilic character of the detergent: the larger the HLB, the more hydrophilic is the detergent. There appears to be some correlation between the HLB value of a detergent and its ability to solubilize membrane proteins. Detergents with a HLB of 12 to 20 are preferred for non-denaturing solubilization of membrane proteins. Detergents with HLBs in the upper end of the above range are preferred for solubilization of extrinsic proteins. Also, the lower the HLB, the more hydrophobic is the detergent and more easily amenable to removal by hydrophobic chromatography.

Summarizing the above properties, it is evident that the performance of a detergent is dependent on the following factors:

| | |
|---------------------------|-------------------------------|
| Detergent concentration | Presence of organic additives |
| Ionic strength | Purity |
| Length of the alkyl chain | Temperature |
| pH | |

Removal of Unbound Detergents

Excess detergent is normally employed in solubilization of membrane proteins. This is to ensure complete dissolution of the membrane and to provide a large number of micelles such that only one protein molecule is present per micelle. However, for further physicochemical and biochemical characterization of membrane proteins, it is often necessary to remove the unbound detergent.

Several methods have been used for detergent removal that take advantage of the general properties of detergents: hydrophobicity, CMC, aggregation number and the charge. The following is a brief description of four commonly used methods.

Hydrophobic Adsorption

This method exploits the ability of detergents to bind to hydrophobic resins. For example, CALBIOCHEM®'s CALBIOSORB™ Adsorbent is a hydrophobic, insoluble resin that can be used in batchwise applications to remove excess detergent. Generally, a solution containing a detergent is mixed with a specific amount of the resin and the mixture is allowed to stand at 4°C or room temperature. The resin with the bound detergent can be removed by centrifugation or filtration. For further details, please refer to Appendix 2. This technique is effective for removal of most detergents. If the adsorption of the protein to the resin

is of concern, the resin can be included in a dialysis buffer and the protein dialyzed. However, this usually requires extended dialyzing periods.

Gel Chromatography

Gel chromatography takes advantage of the difference in size between protein-detergent, detergent-lipid, and homogeneous detergent micelles. In most situations protein-detergent micelles elute in the void volume. The elution buffer should contain a detergent below its CMC value to prevent protein aggregation and precipitation.

Separation by gel chromatography is based on size. Hence, parameters that influence micellar size (ionic strength, pH, and temperature) should be kept constant from experiment to experiment to obtain reproducible results.

Dialysis

When detergent solutions are diluted below the CMC, the micelles are dispersed into monomers. The size of the monomers is usually an order of magnitude smaller than that of the micelles and thus can be easily removed by dialysis. If a large dilution is not practical, micelles can be dispersed by other techniques such as the addition of bile acid salts. For detergents with high CMC, dialysis is usually the preferred choice.

Ion-exchange Chromatography

This method exploits the differences in charge between protein-detergent micelles and protein-free detergent micelles. When non-ionic or zwitterionic detergents are used, conditions can be chosen so that the protein-containing micelles are adsorbed on the ion-exchange resin and the protein-free micelles pass through. Adsorbed protein is washed with detergent-free buffer and is eluted by changing either the ionic strength or the pH. Alternatively, the protein can be eluted with an ionic detergent thus replacing the non-ionic detergent.

Guidelines For Choosing a Detergent

A membrane protein is considered solubilized if it is present in the supernatant after one hour centrifugation of a lysate or a homogenate at 100,000 x g. In most cases, it is also important that the biological activity of the protein be preserved in the supernatant after solubilization by a detergent. Hence, the appropriate detergent is that which yields the maximum amount of biologically active protein in the supernatant.

Given a large number of detergents available today, choosing an appropriate detergent can be a difficult process. Some of the points outlined below can be helpful in selecting a suitable detergent.

1. The first step is a survey of the literature. A detergent that has been used previously for the isolation and characterization of a protein with similar biochemical or enzymological properties should be tried first.
2. Solubility of the detergent at working temperature can be another consideration. For example, ZWITTERGENT® 3-14 is insoluble in water at 4°C while TRITON® X-114 undergoes a phase separation at room temperature.
3. The method of detergent removal can be of an important consideration. If dialysis is to be employed, a detergent with a high CMC is clearly preferred. Alternatively, if ion exchange chromatography is utilized, a non-ionic detergent or a zwittergent is the detergent of choice.
4. Preservation of biological or enzymological activity may require experimenting with several detergents. Not only the type but also the quantity of the detergent used will affect the activity of the protein. For some proteins biological activity is preserved over a very narrow range of concentration of detergent. Below this range the protein is not solubilized and above a particular concentration, the protein is inactivated.
5. Since TRITON® X-100 contains aromatic rings that absorb at 260-280 nm, this detergent should be avoided if the protocols require UV monitoring of protein concentration. Similarly, ionic detergents should be avoided if the proteins are to be separated by isoelectric focusing. For gel filtration of proteins, detergents with smaller aggregation numbers should be considered.
6. Detergents of utmost purity should be used since some detergents such as TRITON® X-100 are generally known to contain peroxides as contaminants. CALBIOCHEM®'s PROTEIN GRADE® or ULTROL® GRADE detergents that have been purified to minimize these oxidizing contaminants should be preferred.

- 7.** CALBIOCHEM® also offers a variety of Molecular Biology Grade detergents for any research where contaminants such as DNase, RNase, and proteases are crucial.
- 8.** A non-toxic detergent should be preferred over a toxic one. For example, digitonin, a cardiac glycoside, should be handled with special care.
- 9.** For as yet unknown reasons, specific detergents often work better for particular isolation procedures. For example, EMPIGEN BB® (Cat. No. 324690) has been found to be the most efficient detergent in solubilization of keratins while preserving their antigenicity. Similarly, *n*-Dodecyl- β -D-maltoside (Cat. No. 324355) has been found to be the detergent of choice for the isolation of cytochrome c oxidase. Hence, some "trial and error" may be required for determining the most optimum conditions for isolation of a membrane protein in its biologically active form.
- 10.** In some cases, it has been observed that the inclusion of non-detergent sulfobetaines (NDSBs) with detergents in the isolation buffer dramatically improves yields of solubilized membrane proteins.

Appendix 1.

Non-Detergent Sulfobetaines

CALBIOCHEM® is pleased to make available a unique line of products, non-detergent sulfobetaines (NDSBs), for protein chemists. NDSBs are zwitterionic compounds. Like zwittergents, NDSBs carry the sulfobetaine hydrophilic head group. However, in contrast to zwittergents, the hydrophobic group in NDSBs is too short for micellar formation even at concentrations as high as 1 M. Hence, they do not behave like detergents. NDSBs were first employed in native isoelectrofocusing gels to neutralize electrostatic interactions without increasing the conductivity. Recently, they have found use in several applications including isolation of membrane proteins, purification of nuclear and halophilic proteins. Presumably, the contribution from the short hydrophobic groups combined with the charge neutralization ability of the sulfobetaine group result in higher yields of membrane proteins. They have also been used in renaturation and refolding of chemically and thermally denatured proteins. It is hypothesized that the hydrophobic group, although short, interacts with the hydrophobic regions of the protein to prevent aggregation during renaturation. They have been used in renaturation of fusion proteins from inclusion bodies.

NDSBs do not interfere with enzymatic assays involving chromogenic substrates bearing nitrophenyl groups and they do not inhibit the activities of β -galactosidase and alkaline phosphatase. In addition, NDSB-195, NDSB-211, and NDSB-221 do not absorb at 280 nm, making them compatible with protein purification procedures in which the protein concentrations are monitored by measuring absorbance at 280 nm.

| Product | *Cat. No. | M.W. |
|----------|-----------|-------|
| NDSB-195 | 480001 | 195.3 |
| NDSB-201 | 480005 | 201.2 |
| NDSB-211 | 480013 | 211.3 |
| NDSB-221 | 480014 | 221.3 |
| NDSB-256 | 480010 | 257.4 |

References:

Chong, Y., and Chen, H. 2000. *Biotechniques* **29**, 1166; Benetti, P.H., et al. 1998. *Protein Expr. Purif.* **13**, 283; Blisnick, T., et al. 1998. *Eur. J. Biochem.* **252**, 537; Vuillard, L., et al. 1998. *Eur. J. Biochem.* **256**, 128; Ochem, A., et al. 1997. *J. Biol. Chem.* **272**, 29919; Goldberg, M.E., et al. 1996. *Folding & Design* **1**, 21; Vuillard, L., et al. 1995. *Anal. Biochem.* **230**, 290; Vuillard, L., et al. 1995. *Biochem. J.* **305**, 337; Vuillard, L., et al. 1994. *FEBS Lett.* **353**, 294.

*You may download product data sheets for these products from our website:

www.calbiochem.com

Appendix 2

CALBIOSORB™ Adsorbent

CALBIOSORB™ Adsorbent (50 ml)
Cat. No. 206550

**CALBIOSORB™ Adsorbent,
Prepacked Column**
(5 ml resin bed + 5 ml buffer reservoir)
Cat. No. 206552

Solubilization of membranes by detergents is essential for their characterization and reconstitution. However, subsequent removal of detergents, particularly the non-ionic detergents with low CMC values, is difficult to achieve.^{1,2} Dialysis, the most common method of detergent removal, usually requires about 200-fold excess of detergent-free buffer with three to four changes over several days. However, it is ineffective for removal of detergents of low CMC values. In addition, prolonged exposure to detergents during dialysis can damage certain membrane proteins.¹ Gel filtration, another common method for detergent removal, is highly effective in the reconstitution of AChR,³ (Ca²⁺ + Mg²⁺)-ATPase,⁴ and lactose transporters.⁵ However, it gives a broader size distribution of vesicles compared to the dialysis method.⁶ Therefore, an expeditious alternative in reconstitution studies is the prior removal of detergents by using a resin capable of effectively binding nondialyzable detergents of low CMC.

CALBIOCHEM® offers an excellent detergent removal product, CALBIOSORB™ Adsorbent. CALBIOSORB™ Adsorbent is a hydrophobic resin that is processed to eliminate unbound organic contaminants, salts, and heavy metal ions and is especially prepared for the removal of detergents from aqueous media. It is supplied in 100 mM Na₂HPO₄, pH 7.0, containing 0.1% sodium azide and can be easily reequilibrated with any other suitable buffer prior to use.

References

1. Jones, O., et al. 1987. In: *Biological Membranes: A Practical Approach* (Findlay, J., and Evens, W., eds.) pp. 139-177. Oxford, IRL Press,
2. Allen, T., et al. 1980. *Biochim. Biophys. Acta* **601**, 328.
3. Mukerjee, P., 1967. *Adv. Colloid. Interface Sci.* **1**, 241.
4. Andersen, J., et al. 1983. *Eur. J. Biochem.* **134**, 205.
5. Furth, A., 1980. *Anal. Biochem.* **109**, 207.
6. Popt, J., and Changeux, J. 1984. *Physiol. Rev.* **64**, 1162.

**Table 1. Detergent Adsorption Capacity
of CALBIOSORB™ Adsorbent**

| <i>Detergent</i> | <i>Mol. Wt.</i> | <i>Type</i> | <i>Adsorption Capacity (mg detergent/ml resin)</i> |
|---|-----------------|--------------|--|
| Cetyltrimethylammonium Bromide (CTAB) | 364.5 | Cationic | 120 |
| CHAPS | 614.9 | Zwitterionic | 110 |
| Cholic Acid, Sodium Salt | 430.6 | Anionic | 73 |
| <i>n</i> -Dodecyl- β -D-maltoside | 510.6 | Non-ionic | 66 |
| <i>n</i> -Hexyl- β -D-glucopyranoside | 264.3 | Non-ionic | 78 |
| Lauryldimethylamine Oxide | 229.4 | Zwitterionic | 66 |
| <i>n</i> -Octyl- β -D-glucopyranoside | 292.4 | Non-ionic | 132 |
| Sodium Dodecyl Sulfate (SDS) | 288.5 | Anionic | 94 |
| <i>n</i> -Tetradecyl- β -D-maltoside | 538.6 | Non-ionic | 161 |
| TRITON® X-100 Detergent | 647.0 (Av.) | Non-ionic | 157 |
| TWEEN® 20, PROTEIN GRADE® Detergent | 1228.0 (Av.) | Non-ionic | 122 |

Detergent absorption capacities were measured by allowing 1.0 g of buffer-free CALBIOSORB™ Adsorbent to equilibrate at room temperature with an excess of detergent (10 ml of 2.0% in H₂O) for 24 hours, then measuring the amount of unabsorbed detergent remaining in the supernatant by gravimetric analysis.

Protocol for Applications Using CALBIOSORB™ Adsorbent, Prepacked Columns

1. Equilibrate the column with 4 to 5 volumes of the sample buffer (e.g., 20 mM sodium phosphate) to remove any sodium azide.
2. Apply the detergent-protein sample to the column.
3. Protein elution from the column may require several column volumes of buffer and can be monitored by UV absorption.

Protocol for Batch Applications Using CALBIOSORB™ Adsorbent

1. Wash CALBIOSORB™ Adsorbent to remove any sodium azide.
2. Calculate the amount of detergent to be removed. For example, 10 ml of 4 mM CHAPS solution contains 24.6 mg of CHAPS.

3. The amount of CALBIOSORB™ Adsorbent required for detergent removal can be determined by inserting the detergent specific adsorption capacity from Table 1 in the following equation:

$$\text{Amount of CALBIOSORB™ Adsorbent} = \frac{\text{Amount of Detergent (mg)}}{\text{Adsorption Capacity (mg/ml)}}$$

(i.e., 24.6 mg CHAPS requires about 0.22 ml of CALBIOSORB™ Adsorbent slurry)

4. Add CALBIOSORB™ Adsorbent directly to the detergent-protein solution. Incubate for 5 minutes at room temperature or for 45 minutes on ice with occasional gentle agitation.
5. Allow the resin to settle. Decant the detergent-free supernatant containing the protein.
6. Dialysis: CALBIOSORB™ Adsorbent may be added directly to a dialysis buffer to facilitate the removal of detergents with low CMC values and to decrease the time required for dialysis when using detergents with higher CMC values. This method is advantageous in that it prevents the adsorption of proteins by the resin.

☞ *A wide variety of application-specific pH and buffer compositions (e.g., HEPES, MOPS, PIPES, Tris, etc.) may be used.*

Storage and Regeneration

Regeneration: Wash with methanol followed by exhaustive washing with water. Reequilibrate with the desired buffer used in the experiment. (NOTE: exhaustive washing is essential to remove methanol from resin). CALBIOSORB™ Adsorbent columns can be used up to ten times before disposal. Regeneration of prepacked CALBIOSORB™ Adsorbent columns is not recommended.

Storage: Wash the resin with a buffer containing 0.1% sodium azide and refrigerate at 4°C.

☞ *When using either the batch or column method, lower ionic strength buffers may decrease the amount of protein absorption by the resin.*

Table 2. Structure and Classification of Detergents



| Detergent class | General structure | Examples |
|------------------|---|---|
| Alkyl glycosides | $R-O-(CH_2)_x-CH_3$ $R-S-(CH_2)_x-CH_3$ | <p>R = glucose x = 8, <i>n</i>-nonyl-β-D-glucopyranoside x = 7, <i>n</i>-octyl-β-D-glucopyranoside x = 6, <i>n</i>-heptyl-β-D-glucopyranoside x = 5, <i>n</i>-hexyl-β-D-glucopyranoside</p> <p>R = maltose x = 11, dodecyl-β-D-maltoside x = 9, decyl-β-D-maltoside</p> <p>R = glucose, x = 7, octyl-β-D-thioglucopyranoside</p> |
| Bile acids |  | <p>x = H, R = O-Na⁺, sodium deoxycholate x = H, R = NHCH₂CH₂SO₃⁻-Na⁺, sodium taurodeoxycholate x = H, R = NHCH₂CO₂⁻-Na⁺, sodium glycodeoxycholate x = OH, R = O-Na⁺, sodium cholate x = OH, R = NHCH₂CH₂SO₃⁻-Na⁺, sodium taurocholate x = OH, R = NHCH₂CO₂⁻-Na⁺, sodium glycocholate</p> |
| Glucamides | $CH_3(CH_2)_x-C(=O)-N(CH_3)-CH_2-C(OH)(H)-C(OH)(H)-C(OH)(H)-C(OH)(H)-OH$  | <p>x = 8, MEGA-10 x = 7, MEGA-9 x = 6, MEGA-8</p> <p>x = H, Deoxy Big CHAP x = OH, Big CHAP</p> |

Table 2. Structure and Classification of Detergents (continued)







| | | |
|--|---|---|
| <p>Polyoxyethylenes, monodisperse and polydisperse</p> |   $\text{CH}_3(\text{CH}_2)_y\text{-O}(\text{CH}_2\text{CH}_2\text{O})_x\text{-H}$ $\text{HO}(\text{CH}_2\text{CH}_2\text{O})_x\text{-(CH}(\text{CH}_3)\text{-CH}_2\text{O)}_y\text{-(CH}_2\text{CH}_2\text{O)}_z\text{-H}$  $W + X + Y + Z = 20$ | <p>x = 9-10, reduced TRITON® X-100 x = 7-8, reduced TRITON® X-114</p> <p>x = 9-10, TRITON® X-100, NP-40 x = 7-8, TRITON® X-114</p> <p>y = 12, X = 8, GENAPOL® X-080 y = 12, X = 10, GENAPOL® X-100 y = 11, x = 8, C₁₂E₈ y = 11, x = 9, C₁₂E₉, THESIT®, LUBROL® PX y = 11, x = 10, GENAPOL® C-100 y = 11, x = 23, BRIJ® 35</p> <p>X = 98, Y = 67, Z = 98, PLURONIC® F-127®</p> <p>R = C₁₁H₂₃CO₂⁻ (laurate), TWEEN® 20 R = C₁₇H₃₃CO₂⁻ (oleate), TWEEN® 80</p> |
|--|---|---|

Table 2. Structure and Classification of Detergents (continued)

Zwittergents

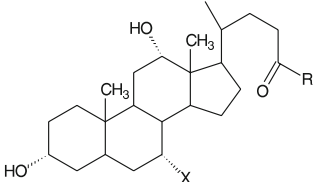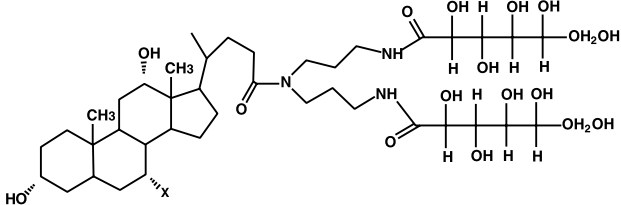

EMPIGEN BB® (*n*-dodecyl-N,N-dimethylglycine)

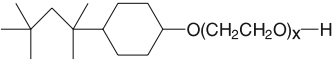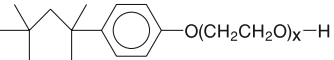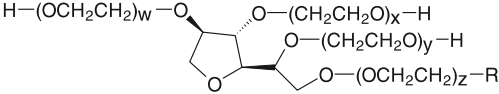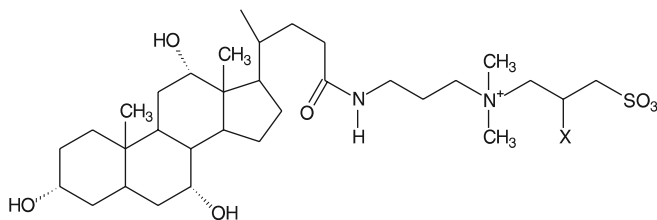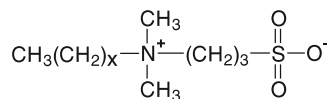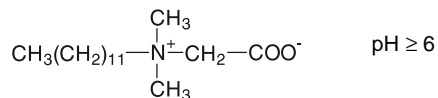x = 7, ZWITTERGENT® 3-08
 x = 9, ZWITTERGENT® 3-10
 x = 11, ZWITTERGENT® 3-12
 x = 13, ZWITTERGENT® 3-14
 x = 15, ZWITTERGENT® 3-16

x = H, CHAPS
 x = OH, CHAPSO

Selected Bibliography

General Properties of Detergents

Amons, R., and Schrier, P.I. 1981. Removal of sodium dodecyl sulfate from proteins and peptides by gel filtration. *Anal. Biochem.* **116**, 439.

Ashani, Y., and Catravas, G.N. 1980. Highly reactive impurities in TRITON® X-100 and BRIJ® 35: partial characterization and removal. *Anal. Biochem.* **109**, 55.

Banerjee, P., et al. 1993. Differential solubilization of membrane lipids by detergents: coenrichment of the sheep brain serotonin 5-HT_{1A} receptor with phospholipids containing predominantly saturated fatty acids. *Arch. Biochem. Biophys.* **305**, 68.

Brito, R.M., and Vaz, W.L.C. 1986. Determination of the critical micelle concentration of surfactants using the fluorescent probe N-phenyl-1-naphthylamine. *Anal. Biochem.* **152**, 250.

Chang, H.W., and Bock, E. 1980. Pitfalls in the use of commercial non-ionic detergents for the solubilization of integral membrane proteins: sulfhydryl oxidizing contaminants and their elimination. *Anal. Biochem.* **104**, 112.

Chattopadhyay, A., and London, E. 1984. Fluorimetric determination of critical micelle concentration avoiding interference from detergent charge. *Anal. Biochem.* **139**, 408.

Furth, A.H., et al. 1984. Separating detergents from proteins. *Methods Enzymol.* **104**, 318.

Helenius, A., and Simons, K. 1975. Solubilization of membranes by detergents. *Biochim. Biophys. Acta* **415**, 29.

Helenius, A., et al. 1979. Properties of detergents. *Methods Enzymol.* **56**, 734.

Hjelmeland, L.M., and Chrambach, A. 1984. Solubilization of functional membrane proteins. *Methods Enzymol.* **104**, 305.

Horigome, T., and Sugano, H. 1983. A rapid method for removal of detergents from protein solutions. *Anal. Biochem.* **130**, 393.

Lever, M. 1977. Peroxides in detergents as interfering factors in biochemical analysis. *Anal. Biochem.* **83**, 274.

Midura, R.J., and Yanagishita, M. 1995. Chaotropic solvents increase the critical micellar concentrations of detergents. *Anal. Biochem.* **228**, 318.

Neugebauer, J. M. 1990. Detergents: An overview. *Methods. Enzymol.* **182**, 239.

Racker, E., et al. 1979. Reconstitution, a way of biochemical research: some new approaches to membrane-bound enzymes. *Arch. Biochem. Biophys.* **198**, 470.

Robinson, N.C., et al. 1984. Phenyl-sepharose mediated detergent exchange chromatography: its application to exchange of detergents bound to membrane proteins. *Biochemistry* **23**, 6121.

Singer, S.J., and Nicolson, G. L. 1972. The fluid mosaic model of the structure of cell membranes. *Science* **175**, 720.

Slinde, E., and Flatmark, T. 1976. Effect of the hydrophile-lipophile balance of non-ionic detergents on the solubilization of biological membranes and their integral b-type cytochromes. *Biochim. Biophys. Acta* **455**, 796.

Storm, D. R., et al. 1976. The HLB dependency for detergent solubilization of hormonally sensitive adenylate cyclase. *J. Supramol. Struct.* **4**, 221.

Tanford, C. 1980. *The Hydrophobic Effect: Formation of Micelles and Biological Membranes* (2nd ed.), New York, Wiley.

Umbreit, J.N., and Strominger, J.L. 1973. Relation of detergent HLB number to solubilization and stabilization of D-alanine carboxypeptidase from *Bacillus subtilis* membranes. *Proc. Natl. Acad. Sci. USA* **70**, 2997.

Detergent-specific References

Zwittergents

Abdullah, K.M., et al. 1995. Purification of baculovirus-overexpressed cytosolic phospholipase A₂ using a single-step affinity column chromatography. *Protein Expr. Purif.* **6**, 291.

Banerjee, P., et al. 1995. Differential solubilization of lipids along with membrane proteins by different classes of detergents. *Chem. Phys. Lipids* **77**, 65.

Cornelius, F., and Skou, J.C. 1984. Reconstitution of sodium-potassium ATPase into phospholipid vesicles with full recovery of its specific activity. *Biochim. Biophys. Acta* **772**, 357.

Fiedler, K., et al. 1993. Glycosphingolipid-enriched, detergent-insoluble complexes in protein sorting in epithelial cells. *Biochemistry* **32**, 6365.

Fricke, B., et al. 2000. Quantitative determination of Zwitterionic detergents using salt-induced phase separation of Triton X-100. *Anal. Biochem.* **281**, 144-50.

Fulop, M.J., et al. 1992. Use of a zwitterionic detergent for the restoration of the antibody binding capacity of immunoblotted *Francisella tularensis* lipopolysaccharide. *Anal. Biochem.* **203**, 141.

Hansel, A., et al. 1994. Isolation and characterization of porin from the outer membrane of *Synechococcus* PCC 6301. *Arch. Microbiol.* **161**, 163.

Hassanain, H.H., et al. 1993. Enhanced gel mobility shift assay for DNA-binding factors. *Anal. Biochem.* **213**, 162.

Lizuka, M., and Fukuda, K. 1993. Purification of the bovine nicotinic acetylcholine receptor α -subunit expressed in baculovirus-infected insect cells. *J. Biochem. (Tokyo)* **114**, 140.

Lowthert, L.A., et al. 1995. Empigen BB: a useful detergent for solubilization and biochemical analysis of keratins. *Biochem. Biophys. Res. Commun.* **206**, 370.

Nguyen, T.D., et al. 1995. Solubilization of receptors for pancreatic polypeptide from rat liver membranes. *Am. J. Physiol.* **268**, G215.

Nollstadt, K.H., et al. 1989. Potential of the sulfobetaine detergent ZWITTERGENT® 3-12 as a desorbing agent in biospecific and bioselective affinity chromatography. *J. Chromatogr.* **497**, 87.

Paik, S.R., et al. 1993. The TF1-ATPase and ATPase activities of assembled alpha 3 beta 3 gamma, alpha 3 beta 3 gamma delta, and alpha 3 beta 3 gamma epsilon complexes are stimulated by low and inhibited by high concentrations of rhodamine 6G whereas the dye only inhibits the alpha 3 beta 3, and alpha 3 beta 3 delta complexes. *J. Bioenerg. Biomembr.* **25**, 679.

Rabilloud, T., et al. 1990. Amidosulfobetaines, a family of detergents with improved solubilization properties: Application for isoelectric focusing under denaturing conditions. *Anal Biochem.* **185**, 94.

Rhinehart-Jones, T., and Greenwalt, D.E. 1996. A detergent-sensitive 113-kDa conformer/complex of CD36 exists on the platelet surface. *Arch. Biochem. Biophys.* **326**, 115.

Riccio, P., et al. 1994. A new detergent to purify CNS myelin basic protein isoforms in lipid-bound form. *NeuroReport* **5**, 689.

Russell-Harde, D., et al. 1995. The use of ZWITTERGENT® 3-14 in the purification of recombinant human β -interferon Ser 17(Betaseron). *J. Interferon Cytokine Res.* **15**, 31.

Schurholz, T., et al. 1992. Functional reconstitution of the nicotinic acetylcholine receptor by CHAPS dialysis depends on the concentrations of salt, lipid, and protein. *Biochemistry* **31**, 5067.

Schurholz, T. 1996. Critical dependence of the solubilization of lipid vesicles by the detergent CHAPS on the lipid composition. Functional reconstitution of the nicotinic acetylcholine receptor into preformed vesicles above the critical micellization concentration. *Biophys. Chem.* **58**, 87.

Spivak, J.L., et al. 1996. Isolation of the full-length murine erythropoietin receptor using a baculovirus expression system. *Blood* **87**, 926.

Stabel, T.J., et al. 1994. Periplasmic location of *Brucella abortus* Cu/Zn superoxide dismutase. *Vet. Microbiol.* **38**, 307.

Stark, R.E., et al. 1984. Physical studies of CHAPS, a new detergent for the study of visual pigments. *J. Phys. Chem.* **88**, 6063.

Valerio, M., and Haraux, F. 1993. Catalytic and activating protons follow different pathways in the H⁺-ATPase of potato tuber mitochondria. *FEBS Lett.* **336**, 83.

Valerio, M., et al. 1993. The electrochemical-proton-gradient-activated states of F₀F₁ ATPase in plant mitochondria as revealed by detergents. *Eur. J. Biochem.* **216**, 565.

Wallace, A.V., and Kuhn, N.J. 1986. Incorporation into phospholipid vesicles of pore-like properties from Golgi membranes of lactating rat mammary gland. *Biochem. J.* **236**, 91.

Warren, B.S., et al. 1996. Purification and stabilization of transcriptionally active glucocorticoid receptor. *J. Biol. Chem.* **271**, 11434.

Xin, H.B., et al. 1995. Affinity purification of the ryanodine receptor/calcium release channel from fast twitch skeletal muscle based on its tight association with FKBP12. *Biochem. Biophys. Res. Commun.* **214**, 263.

Non-ionic Detergents

Bass, W.T., and Bricker, T.M. 1988. Dodecylmaltoside-sodium dodecylsulfate two-dimensional polyacrylamide gel electrophoresis of chloroplast thylakoid membrane proteins. *Anal. Biochem.* **171**, 330.

Begona-Ruiz, M., et al. 1994. An assessment of the biochemical applications of the non-ionic surfactant HECAMEG. *Biochim. Biophys. Acta* **1193**, 301.

Blochet, J.E., et al. 1993. Complete amino acid sequence of puroindoline, a new basic and cystine-rich protein with a unique tryptophan-rich domain, isolated from wheat endosperm by TRITON® X-114 phase partitioning. *FEBS Lett.* **329**, 336.

Dudek, R., et al. 1993. Effect of amphiphiles on nitric oxide synthase in endothelial cells. *Pharmacology* **48**, 374.

Dudeja, P.K., et al. 1995. Reversal of multidrug resistance phenotype by surfactants: relationship to membrane lipid fluidity. *Arch. Biochem. Biophys.* **319**, 309.

El-Kebbi, I.M., et al. 1994. Regulation of the GLUT1 glucose transporter in cultured myocytes: Total number and subcellular distribution as determined by photoaffinity labeling. *Biochem. J.* **301**, 35.

Englund, A.K., et al. 1995. Capillary and rotating-tube isoelectric focusing of a transmembrane protein, the human red cell glucose transporter. *J. Chromatogr. A* **711**, 217.

Florke, R.R., et al. 1993. Differential insertion of insulin receptor complexes into TRITON® X-114 bilayer membranes: evidence for a differential accessibility of membrane-exposed receptor domain. *Eur. J. Biochem.* **211**, 241.

Franek, K.J., et al. 1994. Reliable method for the simultaneous detection of cytoplasmic and surface CD3 epsilon expression by murine lymphoid cells. *Cytometry* **17**, 224.

- Ghebeh, H., et al. 1998. Development of an assay for the measurement of the surfactant Pluronic F-68 in mammalian cell culture medium. *Anal. Biochem.* **262**, 39-44.
- Izawa, S., et al., 1993. Introduction of a series of alkyl thiomaltosides, useful new non-ionic detergents, to membrane biochemistry. *J. Biochem* **113**, 573.
- Kempf, A.C., et al. 1995. Truncated human P450 2D6: expression in *Escherichia coli*, Ni²⁺-chelate affinity purification, and characterization of solubility and aggregation. *Arch. Biochem. Biophys.* **321**, 277.
- Konig, N., and Zampighi, G.A. 1995. Purification of bovine lens cell-to cell channels composed of connexin 44 and connexin 50. *J. Cell. Sci.* **108**, 3091.
- Konrad, R.J., et al. 1995. The heterotrimeric G-protein G_i is localized to the insulin secretory granules of β -cells and is involved in insulin exocytosis. *J. Biol. Chem.* **270**, 12869.
- Lopez-Nicholas, J.M., et al. 1994. An octaethylene glycol monododecyl ether-based mixed micellar assay for lipoxygenase acting at neutral pH. *Anal. Biochem.* **221**, 410.
- Mattsson, J.P., et al. 1994. Isolation and reconstitution of a vacuolar-type proton pump of osteoclast membranes. *J. Biol. Chem.* **269**, 24979.
- Mimura, K., et al. 1993. Change in oligomeric structure of solubilized Na⁺/K⁺-ATPase induced by octaethylene glycol dodecyl ether, phosphatidylserine and ATP. *Biochim. Biophys. Acta* **1145**, 63.
- Moller, J.V., and le Maire, M. 1993. Detergent binding as a measure of hydrophobic surface area of integral membrane proteins. *J. Biol. Chem.* **268**, 18659.
- Nederlof, P.M., et al. 1995. Nuclear localization signals of human and *Thermoplasma proteasomal* subunits are functional *in vitro*. *Proc. Natl. Acad. Sci. USA* **92**, 12060.
- Nguyen, G., and Kruithof, E.K. 1993. A quantitative receptor assay using TRITON[®] X-114 for plasminogen activator binding proteins in solubilized membranes from human liver and placenta. *Anal. Biochem.* **208**, 277.
- Nock, B., et al. 1993. Extracti-Gel D chromatography is a simple, efficient method of removing digitonin during receptor purification: Application to the κ_1 opioid receptor. *J. Neurosci. Methods* **50**, 353.

Ogiso, T., et al. 1994. Mechanism of enhancement effect of *n*-octyl- β -D-thioglycoside on the transdermal penetration of fluorescein isothiocyanate-labeled dextrans and the molecular weight dependence of water-soluble penetrants through stripped skin. *J. Pharm. Sci.* **83**, 1676.

Okamura, S., and Yamashita, S. 1994. Purification and characterization of phosphatidylcholine phospholipase D from pig lung. *J. Biol. Chem.* **269**, 31207.

Pierre, Y., et al. 1995. Purification and characterization of the cytochrome b6f complex from *Chlamydomonas reinhardtii*. *J. Biol. Chem.* **270**, 29342.

Previati, M., et al. 1994. Diacylglycerol kinase activity in rat liver nuclei. *Cell Signal* **6**, 393.

Ramsby, M.L., et al. 1994. Differential detergent fractionation of isolated hepatocytes: biochemical, immunological and two-dimensional gel electrophoresis characterization of cytoskeletal and noncytoskeletal compartments. *Electrophoresis* **15**, 265.

Sivaprasadarao, A., et al. 1994. Solubilization and purification of the retinol-binding protein receptor from human placental membranes. *Biochem. J.* **302**, 245.

Slowiejko, D.M., et al. 1994. Sequestration of muscarinic cholinergic receptors in permeabilized neuroblastoma cells. *J. Neurochem.* **62**, 1795.

Soulimane, T., et al. 1995. Three-dimensional crystals of cytochrome-c oxidase from *Thermus thermophilus* diffracting to 3.8 Å resolution. *FEBS Lett.* **368**, 132.

Strancar, A., et al. 1994. Extraction of TRITON® X-100 and its determination in virus-inactivated human plasma by the solvent-detergent method. *J. Chromatogr. A.* **658**, 475.

Temkin, R.J., et al. 1993. Advantages of digitonin extraction to reveal the intracellular structure of rat glomerular podocytes for high resolution scanning electron microscopy. *Microsc. Res. Tech.* **26**, 260.

Ti, Z.C., et al. 1990. Purification of a membrane glycoprotein with an inositol-containing phospholipid anchor from *Dictyostelium discoideum*. *J. Biotechnol.* **16**, 233.

Virta, M., et al. 1995. Real-time measurement of cell permeabilization with low-molecular weight membranolytic agents. *J. Antimicrob. Chemother.* **36**, 303.

Wallace, P.G., et al. 1994. A method for the determination of the cellular phosphorylation potential and glycolytic intermediates in yeast. *Anal. Biochem.* **222**, 404.

Wong, P. 1993. The state of association of Band 3 of the human erythrocyte membrane: Evidence of a hexamer. *Biochem Biophys Acta* **1151**, 21.

Zardeneta, G., and Horowitz, P.M. 1992. Micelle-assisted protein folding. Denatured rhodanese binding to cardiolipin-containing lauryl maltoside micelles results in slower refolding kinetics but greater enzyme reactivation. *J. Biol. Chem.* **267**, 5811.

Ionic Detergents

Alba, F., et al. 1995. Properties of rat brain dipeptidyl aminopeptidases in the presence of detergents. *Peptides* **16**, 325.

Alba, F., et al. 1995. Comparison of soluble and membrane-bound pyroglutamyl peptidase I activities in rat brain tissues in the presence of detergents. *Neuropeptides* **29**, 103.

Almog, R., et al. 1990. A methodology for determination of Phospholipids. *Anal. Biochem.* **188**, 237.

Bhavsar, J.H., et al. 1994. A method to increase efficiency and minimize anomalous electrophoretic transfer in protein blotting. *Anal. Biochem.* **221**, 234.

Camilleri, P., et al. 1995. High resolution and rapid analysis of branched oligosaccharides by capillary electrophoresis. *Anal. Biochem.* **230**, 115.

Hassaan, A.M., et al. 1995. Calreticulin is the major Ca²⁺ storage protein in the endoplasmic reticulum of the pea plant (*Pisum sativum*). *Biochem. Biophys. Res. Commun.* **211**, 54.

Iwasaki, Y., et al. 1994. Purification and properties of phosphatidylinositol-specific phospholipase C from *Streptomyces antibioticus*. *Biochim. Biophys. Acta.* **1214**, 221.

Kantorow, M., et al. 1995. Conversion from oligomers to tetramers enhances autophosphorylation b lens α -A-crystallin. Specificity between α -A- and α -B-crystallin subunits. *J. Biol. Chem.* **270**, 17215.

Kapp, O.H., and Vinogradov, S.N. 1978. Removal of sodium dodecyl sulfate from Proteins. *Anal. Biochem.* **91**, 230-235.

Komuro, T., et al. 1993. Detection of low molecular size lipopolysaccharide contaminated in dialysates used for hemodialysis therapy with polyacrylamide gel electrophoresis in the presence of sodium deoxycholate. *Int. J. Artif. Organs* **16**, 245.

Muller G., et al. 1994. 4'-Aminobenzamidotaurocholic acid selectively solubilizes glycosyl-phosphatidylinositol-anchored membrane proteins and improves lipolytic cleavage of their membrane anchors by specific phospholipases. *Arch. Biochem. Biophys.* **309**, 329.

Palmer, M., et al. 1995. Kinetics of streptolysin O self-assembly. *Eur. J. Biochem.* **231**, 388.

Rozema, Z., and Gellman, S.H. 1996. Artificial chaperone-assisted refolding of carbonic anhydrase. *J. Biol. Chem.* **271**, 3478.

Siler, D.J., and Cornish, K. 1995. Measurement of protein in natural rubber latex. *Anal. Biochem.* **229**, 278.

Spivak, W., et al. 1988. Spectrophotometric determination of the critical micellar concentration of bile salts using bilirubin monoglucuronide as a micellar probe. Utility of derivative spectroscopy. *Biochem. J.* **252**, 275.

Sundquist, B., et al. 1983. Assay of detergents by rocker electrophoresis in agarose gels containing red blood cells: "Rocker hemolysis". *Biochem. Biophys. Res. Commun.* **114**, 699.

Tadey, T., and Purdy, W.C. 1993. Effect of detergents on the electrophoretic behavior of plasma apolipoproteins in capillary electrophoresis. *J. Chromatogr. A* **652**, 131.

Taipale, J., et al. 1995. Human mast cell chymase and leukocyte elastase release latent transforming growth factor- β 1 from the extracellular matrix of cultured human epithelial and endothelial cells. *J. Biol. Chem.* **270**, 4689.

Non-ionic Detergents

| Product | Cat. No. | M.W. (anhydrous) | CMC ^b (mM) | Aggregation No. | Average Micellar Weight | size |
|---|----------|---------------------|--------------------------|-----------------|----------------------------|----------------------|
| APO-10 | 178375 | 218.3 | 4.6 | 131 | 28,000 | 1 g |
| APO-12 | 178377 | 246.4 | 0.568 | 2232 | 500,000 | 1 g |
| Big CHAP | 200965 | 878.1 | 3.4 | 10 | 8,800 | 1 g |
| Big CHAP, Deoxy | 256455 | 862.1 | 1.1 - 1.4 | 8-16 | 10,500 | 250 mg 1 g |
| BRIJ® 35, PROTEIN GRADE® Detergent, 30% Solution | 203724 | — | 0.09 | 40 | 48,000 | 100 ml 1 L |
| BRIJ® 35, PROTEIN GRADE® Detergent, 10% Solution, Sterile-Filtered | 203728 | — | 0.09 | 40 | 48,000 | 50 ml |
| C ₁₂ E ₅ | 205524 | — | — | — | — | 1 g |
| C ₁₂ E ₆ | 205527 | 450.7 | 0.087 | — | — | 1 g |
| C ₁₂ E ₈ | 205528 | 538.8 | 0.11 | 123 | 66,000 | 1 g |
| C ₁₂ E ₉ | 205529 | 582.8 | 0.08 | — | — | 1 g |
| Cyclohexyl- <i>n</i> -ethyl-β-D-maltoside, ULTROL® Grade | 239774 | 452.5 | 120 | — | — | 1 g |
| Cyclohexyl- <i>n</i> -hexyl-β-D-maltoside, ULTROL® Grade | 239775 | 508.6 | 0.56 | 63 | 32,000 | 500 mg 1 g |
| Cyclohexyl- <i>n</i> -methyl-β-D-maltoside, ULTROL® Grade | 239776 | 438.5 | 340 | — | — | 1 g |
| <i>n</i> -Decanoylsucrose | 252721 | 496.6 | 2.5 | — | — | 1 g 5 g |
| <i>n</i> -Decyl-β-D-maltopyranoside, ULTROL® Grade | 252718 | 482.6 | 1.6 | — | — | 1 g 5 g |
| <i>n</i> -Decyl-β-D-thiomaltoside, ULTROL® Grade | 252725 | 498.6 | 0.9 | — | — | 500 mg |
| Digitonin, High Purity | 300410 | 1229.3 | — | 5 - 6 | 7000 | 250 mg 1 g 5 g |

Non-ionic Detergents cont.

| Product | Cat. No. | M.W. (anhydrous) | CMC ^b (mM) | Aggregation No. | Average Micellar Weight | size |
|--|----------|---------------------|--------------------------|-----------------|----------------------------|------------------------------|
| Digitonin, Alcohol-Soluble, High Purity | 300411 | 1229.3 | — | 5 - 6 | 7000 | 250 mg 1 g |
| <i>n</i> -Dodecanoylsucrose | 324374 | 524.6 | 0.3 | — | — | 1 g 5 g |
| <i>n</i> -Dodecyl- β -D-glucopyranoside | 324351 | 348.5 | 0.13 | 200 | 70,000 | 1 g |
| <i>n</i> -Dodecyl- β -D-maltoside, ULTROL [®] Grade | 324355 | 510.6 | 0.1-0.6 | 98 | 50,000 | 500 mg 1 g 5 g 25 g |
| ELUGENT™ Detergent, 50% Solution | 324707 | — | — | — | — | 100 ml |
| GENAPOL [®] C-100, PROTEIN GRADE [®] Detergent, 10% Solution | 345794 | 627.0 | — | — | — | 50 ml |
| GENAPOL [®] X-80, PROTEIN GRADE [®] Detergent, 10% Solution | 345796 | 553.0 | 0.06 - 0.15 | — | — | 50 ml |
| GENAPOL [®] X-100, PROTEIN GRADE [®] Detergent, 10% Solution | 345798 | 641.0 | 0.15 | 88 | 56,000 | 50 ml |
| <i>n</i> -Heptyl- β -D-glucopyranoside | 375655 | 278.3 | 79 | — | — | 1 g |
| <i>n</i> -Heptyl- β -D-thioglucopyranoside, ULTROL [®] Grade, 10% Solution | 375659 | 294.4 | 30 | — | — | 10 ml 50 ml |
| <i>n</i> -Hexyl- β -D-glucopyranoside | 376965 | 264.3 | 250 | — | — | 1 g |
| MEGA-8, ULTROL [®] Grade | 444926 | 321.5 | 58 | — | — | 1 g 5 g |
| MEGA-9, ULTROL [®] Grade | 444930 | 335.5 | 19 - 25 | — | — | 1 g 5 g |
| MEGA-10, ULTROL [®] Grade | 444934 | 349.5 | 6 - 7 | — | — | 1 g 5 g |

Non-ionic Detergents cont.

| Product | Cat. No. | M.W. (anhydrous) | CMC ^b (mM) | Aggregation No. | Average Micellar Weight | size |
|---|----------|---------------------|--------------------------|-----------------|----------------------------|-----------------------------|
| <i>n</i> -Nonyl- β -D-glucopyranoside | 488285 | 306.4 | 6.5 | — | — | 1 g |
| NP-40 | 492015 | 603.0 | 0.05-0.3 | — | — | 100 ml 500 ml 1000 ml |
| NP-40, PROTEIN GRADE [®] Detergent, 10% Solution | 492017 | 603.0 | 0.05-0.3 | — | — | 50 ml |
| <i>n</i> -Octanoyl- β -D-glucosylamine (NOGA) | 488100 | 305.4 | 80 | — | — | 500 mg 1 g 5 g |
| <i>n</i> -Octanoylsucrose | 494466 | 468.5 | 24.4 | — | — | 1 g 5 g |
| <i>n</i> -Octyl- β -D-glucopyranoside | 494459 | 292.4 | 20 - 25 | 84 | 25,000 | 1 g 5 g 25 g |
| <i>n</i> -Octyl- β -D-glucopyranoside, ULTROL [®] Grade | 494460 | 292.4 | 20 - 25 | 84 | 25,000 | 250 mg 1 g 5 g |
| <i>n</i> -Octyl- β -D-maltopyranoside | 494465 | 454.5 | 23.4 | 84 | 38,000 | 1 g |
| <i>n</i> -Octyl- β -D-thioglycopyranoside, ULTROL [®] Grade | 494461 | 308.4 | 9 | — | — | 5 g 25 g |
| PLURONIC [®] F-127, PROTEIN GRADE [®] Detergent, 10% Solution | 540025 | — | 4 - 11 | — | — | 50 ml |
| TRITON [®] X-100 | 648462 | 625 (avg.) | 0.2 - 0.9 | 100 - 155 | 80,000 | 1 kg |
| TRITON [®] X-100, PROTEIN GRADE [®] Detergent, 10% Solution | 648463 | 625 (avg.) | 0.2 - 0.9 | 100 - 155 | 80,000 | 50 ml |
| TRITON [®] X-100, Molecular Biology Grade | 648466 | 625 (avg.) | 0.2 - 0.9 | 100 - 155 | 80,000 | 50 ml |
| TRITON [®] X-100, Hydrogenated | 648465 | 631 (avg.) | 0.25 | 100 - 155 | 80,000 | 10 g |

Non-ionic Detergents cont.

| Product | Cat. No. | M.W. (anhydrous) | CMC ^b (mM) | Aggregation No. | Average Micellar Weight | size |
|---|----------|---------------------|--------------------------|-----------------|----------------------------|--------|
| TRITON® X-100, Hydrogenated, PROTEIN GRADE® Detergent, 10% Solution | 648464 | 631 (avg.) | 0.25 | 100 - 155 | 80,000 | 10 ml |
| TRITON® X-114, PROTEIN GRADE® Detergent, 10% Solution | 648468 | 537 (avg.) | 0.35 | — | — | 50 ml |
| TWEEN® 20 | 655205 | 1228 (avg.) | 0.059 | — | — | 250 ml |
| TWEEN® 20, Molecular Biology Grade | 655204 | — | 0.059 | — | — | 100 ml |
| TWEEN® 20, PROTEIN GRADE® Detergent, 10% Solution | 655206 | 1228 (avg.) | 0.059 | — | — | 50 ml |
| TWEEN® 80, PROTEIN GRADE® Detergent, 10% Solution | 655207 | 1310 (avg.) | 0.012 | 58 | 76,000 | 50 ml |
| <i>n</i> -Undecyl-β-D-maltoside, ULTROL® Grade | 662085 | 496.6 | 0.59 | — | — | 500 mg |

Ionic Detergents

| Product | Cat. No. | M.W. (anhydrous) | CMC ^b (mM) | Aggregation No. | Average Micellar Weight | size |
|--|----------|---------------------|--------------------------|-----------------|----------------------------|-----------------------|
| BATC | 196950 | 655.8 | — | — | — | 500 mg |
| Cetyltrimethylammonium Bromide (CTAB) Molecular Biology Grade | 219374 | 364.5 | 1.0 | 170 | 62,000 | 100 g |
| Chenodeoxycholic Acid, Free Acid | 2204 | 392.6 | — | — | — | 5 g |
| Chenodeoxycholic Acid, Sodium Salt | 220411 | 414.6 | — | — | — | 1 g 5 g |
| Cholic Acid, Sodium Salt | 229101 | 430.6 | 9 - 15 | 2.0 | 900 | 50 g 250 g 1 kg |
| Cholic Acid, Sodium Salt, ULTROL® Grade | 229102 | 430.6 | 9 - 15 | 2.0 | 900 | 1 g 5 g |
| Deoxycholic Acid, Sodium Salt | 264101 | 414.6 | 4 - 8 | 4-10 | 1600-4100 | 25 g 100 g |
| Deoxycholic Acid, Sodium Salt, ULTROL® Grade | 264103 | 414.6 | 2 - 6 | 3-12 | 1200-4900 | 5 g 25 g 100 g |
| 7 α , 12 α -Dihydroxy-5 β -cholic Acid | 305705 | 392.6 | — | — | — | 100 mg |
| Glycholic Acid, Sodium Salt | 360512 | 487.6 | 7.1 | 2.1 | 1000 | 1 g 5 g |

Ionic Detergents

| Product | Cat. No. | M.W. (anhydrous) | CMC ^b (mM) | Aggregation No. | Average Micellar Weight | size |
|--|----------|---------------------|--------------------------|-----------------|----------------------------|---------------|
| Glycodeoxycholic Acid, Sodium Salt | 361311 | 471.6 | 2.1 | 2.1 | 900 | 5 g |
| Lauroylsarcosine, Sodium Salt | 428010 | 293.4 | — | 2.0 | 600 | 5 g |
| Sodium <i>n</i> -Dodecyl Sulfate (SDS) | 428015 | 288.5 | 7 - 10 | 62 | 18,000 | 1 kg |
| Sodium <i>n</i> -Dodecyl Sulfate (SDS), High Purity | 428016 | 288.5 | 7 - 10 | 62 | 18,000 | 25 g |
| Sodium <i>n</i> -Dodecyl Sulfate (SDS), Molecular Biology Grade | 428023 | 288.5 | 7 - 10 | 62 | 18,000 | 50 g 500 g |
| Sodium <i>n</i> -Dodecyl Sulfate (SDS), 30% Solution | 428018 | 288.5 | 7 - 10 | 62 | 18,000 | 200 ml |
| Taurochenodeoxycholic Acid, Sodium Salt | 580211 | 521.7 | — | — | — | 1 g 5 g |
| Taurocholic Acid, Sodium Salt | 580217 | 537.7 | 3 - 11 | 4 | 2100 | 5 g 25 g |
| Taurocholic Acid, Sodium Salt, ULTROL [®] Grade | 580218 | 537.7 | 3 - 11 | 4 | 2100 | 1 g 5 g |
| Taurodehydrocholic Acid, Sodium Salt | 580219 | 531.6 | — | — | — | 1 g |
| Taurodeoxycholic Acid, Sodium Salt | 580221 | 521.7 | 1 - 4 | 6 | 3100 | 5 g 50 g |
| Tauroolithocholic Acid, Sodium Salt | 580352 | 505.7 | — | — | — | 1 g 5 g |
| Tauroursodeoxycholic Acid, Sodium Salt | 580549 | 521.7 | — | — | — | 1 g 5 g |
| TOPPS | 615000 | 350.5 | 4.5 | — | — | 5 g |

Zwitterionic Detergents

| Product | Cat. No. | M.W. (anhydrous) | CMC ^b (mM) | Aggregation No. | Average Micellar Weight | size |
|--|----------|---------------------|--------------------------|-----------------|----------------------------|----------------------------|
| ASB-14 | 182750 | 434.7 | — | — | — | 5g, 25 g |
| ASB-16 | 182755 | 462.7 | — | — | — | 5g, 25 g |
| CHAPS | 220201 | 614.9 | 6-10 | 4-14 | 6000 | 1 g 5 g 10 g 25 g |
| CHAPSO | 220202 | 630.9 | 8 | 11 | 7000 | 1 g 5 g |
| DDMAB | 252000 | 299.5 | 4.3 | — | — | 5 g |
| DDMAU | 252005 | 397.7 | 0.13 | — | — | 5 g |
| EMPIGEN BB® Detergent, 30% Solution | 324690 | 272.0 | 1.6-2.1 | — | — | 100 ml |
| Lauryldimethylamine Oxide (LDAO), 30% Solution | 428011 | 229.4 | 1-2 | 76 | 17,000 | 100 ml |
| ZWITTERGENT® 3-08 Detergent | 693019 | 279.6 | 330 | — | — | 5 g |
| ZWITTERGENT® 3-10 Detergent | 693021 | 307.6 | 25-40 | 41 | 12,500 | 5 g 25 g |
| ZWITTERGENT® 3-12 Detergent | 693015 | 335.6 | 2-4 | 55 | 18,500 | 5 g 25 g |
| ZWITTERGENT® 3-14 Detergent | 693017 | 363.6 | 0.1-0.4 | 83 | 30,000 | 5 g 25 g 100 g |
| ZWITTERGENT® 3-16 Detergent | 693023 | 391.6 | 0.01-0.06 | 155 | 60,000 | 5 g 25 g |

a. Average molecular weights are given for detergents composed of mixtures of chain lengths; b. Temperature: 20 - 25°C

TRADEMARKS

BRIJ® and **TWEEN®** are registered trademarks of ICI Americas, Inc.

EMPIGEN BB® is a registered trademark of Allbright & Wilson.

GENAPOL® is a registered trademark of Hoechst AG.

ULTROL®, **PROTEIN GRADE®** and **ZWITTERGENT®** are registered trademarks of Calbiochem-Novabiochem Corporation.

CALBIOSORB™ Adsorbent and **ELUGENT™ Detergent** are trademarks of Calbiochem-Novabiochem Corporation

LUBROL® is a registered trademark of Imperial Chemical Inc.

PLURONIC® is a registered trademark of Wyandotte Chemicals Corporation.

TRITON X® is a registered trademark of Rohm and Hass.

*For Fast Efficient Service
Call or Write...*



USA, CANADA, & MEXICO
CALBIOCHEM-NOVABIOCHEM
P.O. Box 12087
La Jolla, CA 92039-2087

To Place Your Order:
800-854-3417
E-mail: orders@calbiochem.com

For Technical Inquiries:
800-628-8470
E-mail: technical@calbiochem.com

Street Address:
10394 Pacific Center Court
San Diego, California
92121
Website: <http://www.calbiochem.com>

You may also order
CALBIOCHEM®
products through:





www.vwrsp.com



A VWR COMPANY

www.vwrcanlab.com

Customer Service
(800) 932-5000

 CALBIOCHEM® A brand of  biosciences, inc.
An Affiliate of Merck KGaA, Darmstadt, Germany

Doc No. **CB0068-401**

GERMANY
CALBIOCHEM-NOVABIOCHEM GmbH
Telephone: +49 (0)6196 - 63955
Freecall: (0800) - 6931 000
Fax: +49 (0)6196 - 62361
Freefax: (0800) 62361 00

E-mail: customer.service@calbiochem-novabiochem.de

MERCK EUROLAB GmbH
Tel: 01805 70 20 00
Fax: 01805 70 22 22
E-mail: biotech@merckeurolab.de

UNITED KINGDOM
CN BIOSCIENCES (UK) LTD.
Freefone: 0800 622935
Telephone: 0115 943 0840
Fax: 0115 943 0951
E-mail:
customer.service@cruk.co.uk
Website: <http://www.cnbiosciences.co.uk>

From Republic of Ireland
Telephone: 1800 409445 Toll free
+44 115 943 0840
Fax: +44 115 943 0951

Merck Eurolab Ltd
Freefone: 0800 223344
Fax: 01455 558586
E-mail: uk.sales@merckeurolab.ltd.uk